



Check pesticide residues in fruits and vegetables with Ringbio Rapid Test Kit

- 🌸 One Step Universal Preparation
- 🌸 Rapid 10min Test
- 🌸 No Special Instrument Required



RING BIOTECHNOLOGY CO LTD



Scan to read PDF copy

Pesticides Residue Rapid Test Kit


Item #	Product	Item #	Product
V001	Carbendazim Rapid Test Kit	V011	Albendazole Rapid Test Kit
V002	Paraquat Rapid Test Kit	V013	Fipronil Rapid Test Kit
V003	Chlorothalonil Rapid Test Kit	V017	Carbendazim & Acetamiprid & Carbaryl 3in1 Rapid Test Kit
V004	Acetamiprid Rapid Test Kit	V018	Thiamethoxam & cyromazine, TC Combo Test Kit
V005	Imidacloprid Rapid Test Kit	V029	Glyphosate test kit
V006	Cyromazine Rapid Test Kit		
V007	Carbaryl Rapid Test Kit		
V008	Thiamethoxam Rapid Test Kit		
V009	Triazolone-Triadimenol Test Kit		
V010	Carbofuran Rapid Test Kit		

** all kits are 48 tests or 96 tests per box. Sample buffer and test strips are included.


Procedures for testing



Take sample piece into a bag and mix with sample buffer

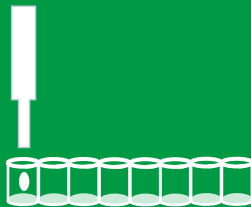

Add 200ul mixed sample into vial mix and incubate for 3min

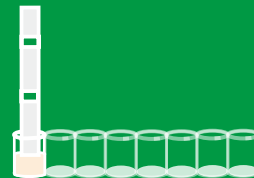

Insert test strip Wait for 5min



Check the result

** all samples shall be homogenized before testing. For colored sample, filtration or centrifuge may require. Please always check the kit instruction manual for detailed operation guide.



Ring Biotechnology Co Ltd

Add: Building 3, Zhongtongtai TechnoPark, BDA, Beijing 100176, China

Tel: +86-10-56267496 / +86-10-26242030

Email: info@ringbio.com

Web: www.ringbio.com